


CIPAC

COLLABORATIVE INTERNATIONAL PESTICIDES ANALYTICAL COUNCIL LIMITED

Commission Internationale des Méthodes d'Analyse des Pesticides (CIMAP)

CHAIRMAN:
Dr Markus.D. Müller
Swiss Federal Research Station
CH-8820 Wädenswil, Switzerland
Telephone: +41 1 783 64 12
Telefax: +41 1 783 64 39
E-Mail: Markus.Mueller@faw.admin.ch

TREASURER:
Mr Brian E. Hocken
17 Claygate Avenue
Harpenden, Herts, AL5 2HE, England
Telephone: +44 1582 764277
Telefax: +44 1582 461625
E-Mail: brian@hockene.freeserve.co.uk

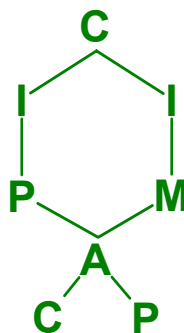

SECRETARY:
Mr László Bura
Central Service for Plant Protection & Soil Conservation
Budaörsi út 141-145
1118 Budapest, Hungary
Telephone: +36 1 309 10 27
Telefax: +36 1 246 29 56
E-Mail: bura.laszlo@ntksz.ontsz.hu

ASSISTANT SECRETARY:
Dr Anna Kouppari
Analytical Laboratories, Dept of Agriculture
Nicosia 1412, Cyprus
Telephone: +357 22 301 889
Telefax: +357 22 781425
E-mail: anna.kouppari@cytanet.com.cy

Internet: <http://www.cipac.org>

CIPAC INFORMATION SHEET No. 241

Your Ref:
Our Ref: LB
29 July 2003

743. prallethrin

Dear Member,

This is to inform you that JAPAC is planning to organize a collaborative trial on a method for prallethrin.

Principle. The content of prallethrin is determined by capillary GC (column 0.25mm ID x 30m, coated with a 0.25µm film of crosslinked 14% cyanopropylphenyl 86% dimethyl polysiloxane e.g. DB-1701) using split injection, flame ionisation detector and triphenyl phosphate as internal standard.

Scope. It is intended to analyze five samples consisting of two technicals and three liquid vaporizers. Prallethrin analytical standard will be supplied.

Time schedule. The collaborative trial is intended to be carried out during October 2003 - January 2004. Those who wish to participate should contact the organizer by 31 August 2003. Participants who perform the trial together with 356 *d*-phenothrin would be welcomed.

Contact:

Dr Ritsuko Furuta
Sumitomo Chemical Co., Ltd.
Organic Synthesis Research Laboratory
Analytical Science Group
3-1-98 Kasugadenaka Konohanaku,
Osaka, 554-8558 Japan

TEL: +81-6-6466-5171
FAX: +81-6-6466-5459
E-mail: furutar@sc.sumitomo-chem.co.jp

Yours sincerely

L. Bura
Secretary