

CIPAC

COLLABORATIVE INTERNATIONAL PESTICIDES ANALYTICAL COUNCIL LIMITED


Commission Internationale des Méthodes d'Analyse des Pesticides (CIMAP)

CHAIRMAN:

Dr. Ralf Hänel
Federal Office of Consumer Protection and Food Safety
Messeweg 11/12
38104 Braunschweig, Germany
Telephone: +49 531 299 3506
Telefax: +49 531 299 3002
E-Mail: ralf.haenel@bvl.bund.de

TREASURER:

Brian E. Hocken
17 Claygate Avenue
Harpenden, Herts, AL5 2HE, England
Telephone: +44 1582 764277
Telefax: +44 1582 461625
E-Mail: brian@hockene.freereserve.co.uk

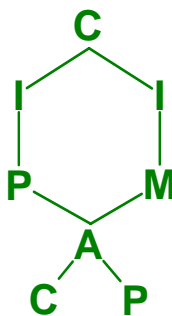

SECRETARY:

László Bura
European Food Safety Authority
Largo N. Palli 5/A, I-43100 Parma, Italy
Telephone: +39 0521 036827
Telefax: +39 0521 0360827
E-Mail: laszlo.bura@efsa.europa.eu

ASSISTANT SECRETARY:

Dr. Eric Sandmann
Division Pesticide Science
Plant Protection Research Institute
Agricultural Research Council
Pretoria, South Africa
Tel +27 12 8088000
Fax +27 12 8088299
E-Mail: SandmannE@arc.agric.za

Internet: <http://www.cipac.org>

CIPAC INFORMATION SHEET No 281

Your Ref:
Our Ref: LB
12 November 2008

738. clothianidin

This is to inform you that JAPAC is planning to organize a collaborative trial on a method for clothianidin.

Principle. The content of clothianidin is determined by reversed phase high performance liquid chromatography on a 150 x 4.6 mm ODS column (Zorbax Eclipse XDB-C18 (5 µm), or equivalent) with mobile phase composed of acetonitrile/water/phosphoric acid, 150/850/1 (v/v) and UV detection at 269 nm. Quantification is done by external standardization.

Scope. It is intended to analyze six samples consisting of one technical material, one GR, one SC, one WG and two SG formulations. Clothianidin analytical standard will be supplied.

Time schedule. The collaborative trial is intended to be carried out during December 2008 - March 2009.

Those who wish to participate should contact the organizer by 30th November 2008.

Contact:

Mr. Yasushi Asada
Sumitomo Chemical Co., Ltd.
Organic Synthesis Research Laboratory
Analytical Science Group
3-1-98 Kasugadenaka Konohanaku,
Osaka, 554-8558 Japan

TEL: +81-6-6466-5171
FAX: +81-6-6466-5459
E-mail: asada@sc.sumitomo-chem.co.jp

Yours sincerely

L. Bura
Secretary